


DEQAS April 2017 - 25-OHD Method Means (+/-1SD) for Major Method Groups\*

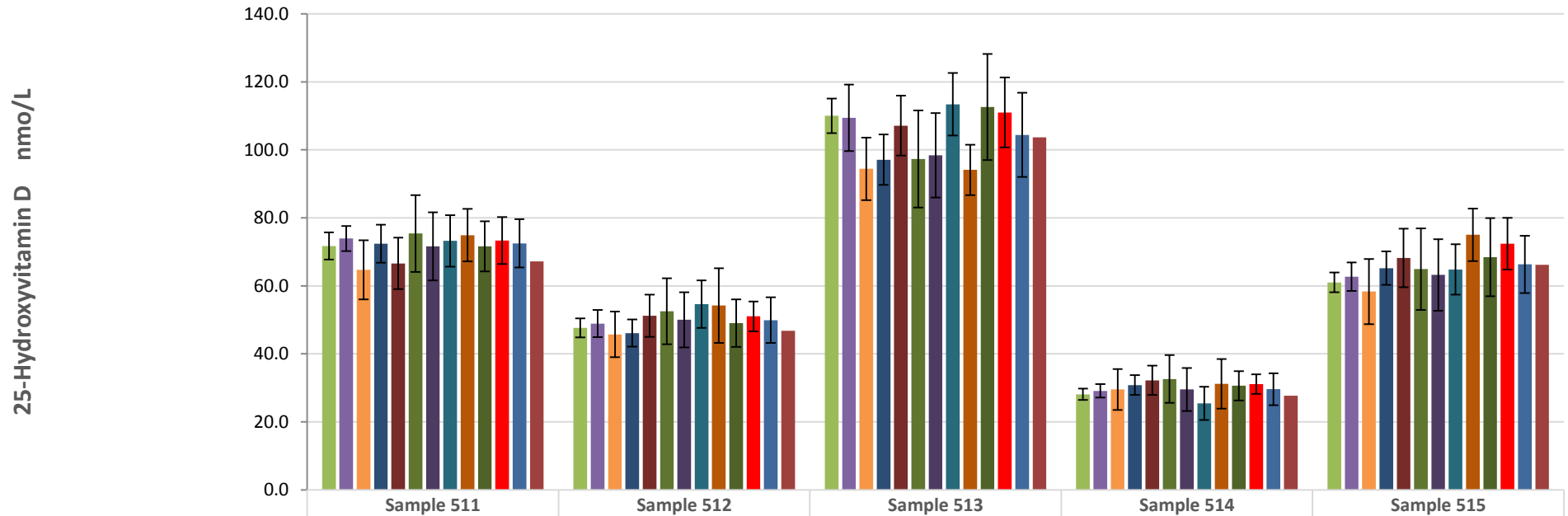

	Sample 511	Sample 512	Sample 513	Sample 514	Sample 515
■ Abbott Architect New Kit (n=67)	71.7	47.6	110.0	28.1	61.0
■ Abbott Architect Old Kit (n=12)	73.9	48.9	109.4	29.1	62.7
■ Beckman Unicel (n=33)	64.7	45.7	94.4	29.5	58.3
■ DiaSorin Liaison (n=211)	72.4	46.1	97.1	30.8	65.2
■ Euroimmun ELISA (n=11)	66.6	51.2	107.1	32.2	68.2
■ IDS EIA (n=12)	75.4	52.5	97.3	32.6	64.9
■ IDS iSYS (n=74)	71.6	50.0	98.4	29.5	63.2
■ Roche Total (n=168)	73.2	54.6	113.4	25.4	64.8
■ Siemens (n=61)	74.9	54.2	94.1	31.2	75.0
■ HPLC (n=18)	71.6	49.0	112.6	30.6	68.4
■ LC-MS/MS (n=159)	73.3	51.0	111.0	31.1	72.4
■ ALTM (n=877)	72.5	49.9	104.4	29.6	66.3
■ TARGET VALUE	67.2	46.8	103.7	27.7	66.2

\* Methods with 10 or more results returned