


DEQAS April 2018 - 25-OHD Method Means (+/-1SD) for Major Method Groups\*

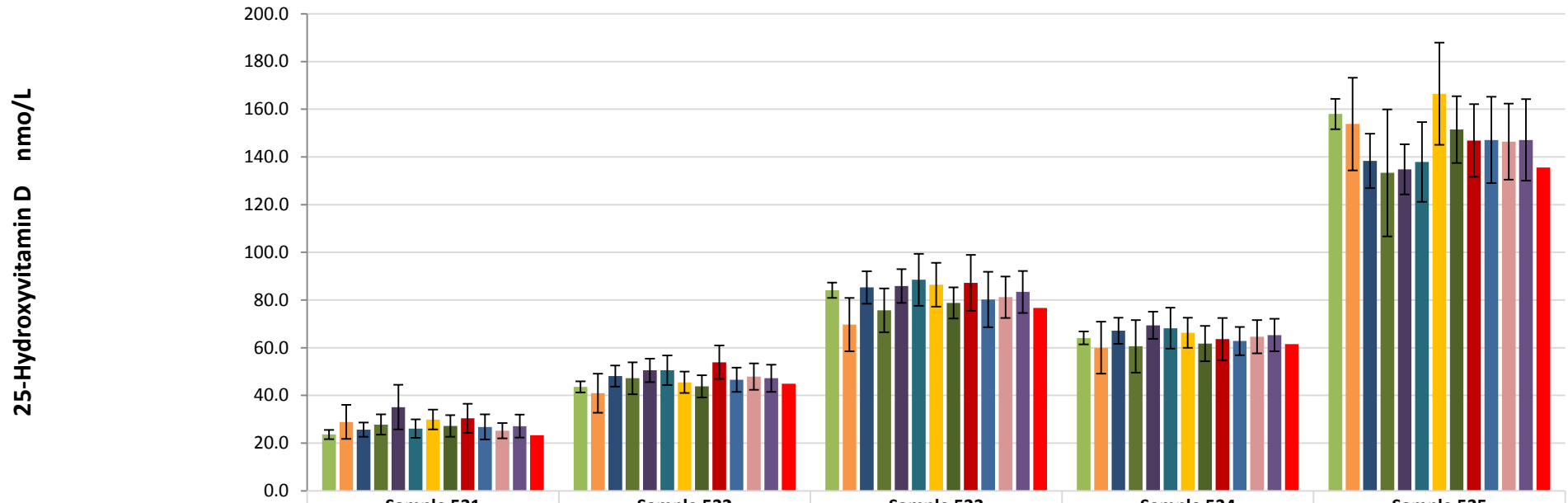

	Sample 531	Sample 532	Sample 533	Sample 534	Sample 535
Abbott Architect New Kit (n=68)	23.6	43.6	84.1	64.1	158.0
Beckman Unicel (n=31)	28.9	40.9	69.7	60.0	153.8
DiaSorin Liaison (n=198)	25.6	48.1	85.3	67.1	138.3
IDS EIA (n=10)	27.8	47.2	75.7	60.6	133.3
IDS iSYS (n=43)	35.1	50.5	85.9	69.4	134.8
IDS-iSYS New (n=16)	26.1	50.5	88.5	68.2	137.9
Roche Total (n=108)	29.9	45.5	86.4	66.3	166.5
Roche Vitamin D total II (n=35)	27.2	43.8	78.8	61.7	151.5
Siemens (n=55)	30.4	53.9	87.2	63.6	146.9
HPLC (n=14)	26.8	46.6	80.2	62.8	147.1
LC-MS/MS (n=142)	25.2	47.9	81.2	64.6	146.4
ALTM (n=776)	27.1	47.2	83.4	65.3	147.1
TARGET VALUE	23.3	44.9	76.7	61.5	135.6

\* Methods with 10 or more results returned