


DEQAS July 2017 - 25-OHD Method Means (+/-1SD) for Major Method Groups*

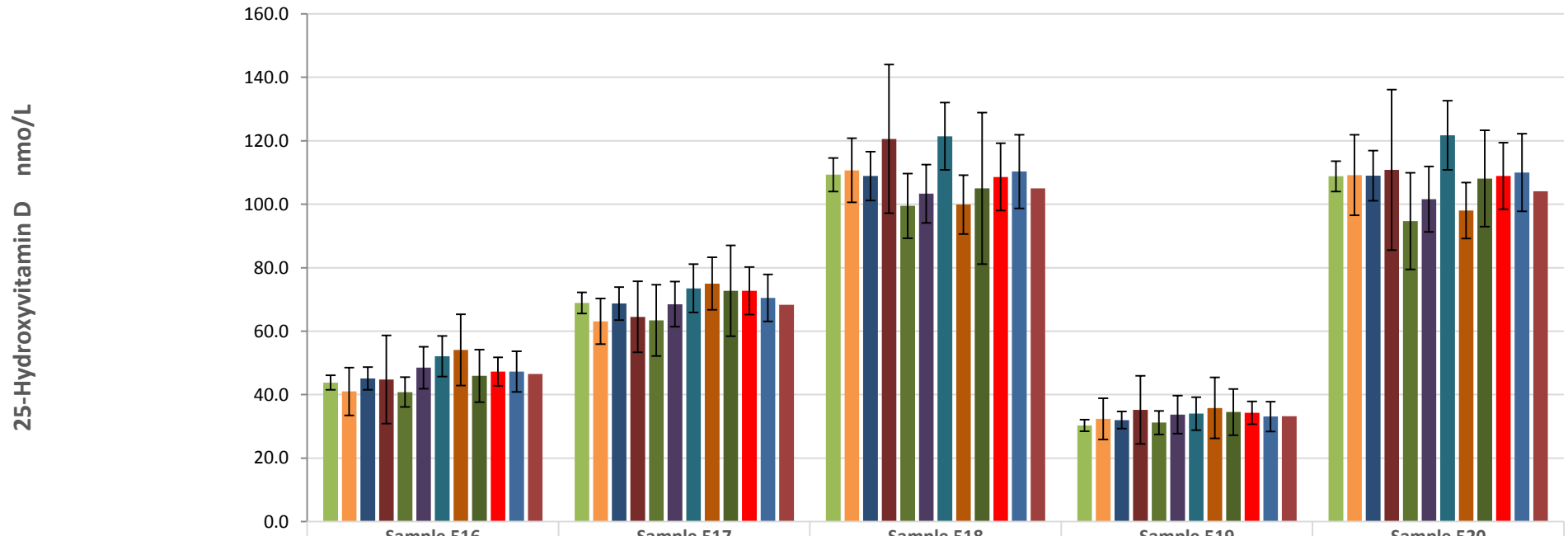

	Sample 516	Sample 517	Sample 518	Sample 519	Sample 520
Abbott Architect New Kit (n=71)	43.8	68.9	109.3	30.3	108.8
Beckman Unicel (n=32)	41.0	63.1	110.7	32.4	109.2
DiaSorin Liaison (n=214)	45.1	68.7	108.9	32.0	109.0
Euroimmun ELISA (n=10)	44.8	64.5	120.6	35.2	110.8
IDS EIA (n=12)	40.8	63.4	99.5	31.2	94.7
IDS iSYS (n=63)	48.5	68.5	103.3	33.7	101.6
Roche Total (n=165)	52.1	73.5	121.4	34.0	121.7
Siemens (n=60)	54.1	75.0	99.9	35.8	98.0
HPLC (n=18)	45.9	72.7	105.0	34.5	108.1
LC-MS/MS (n=159)	47.3	72.7	108.6	34.3	108.9
ALTM (n=867)	47.3	70.5	110.3	33.1	110.0
TARGET VALUE	46.5	68.3	105.0	33.2	104.1

* Methods with 10 or more results returned